

Ordering Details

purple surgical


purple surgical

Ultimate Endoscopic Linear Cutter Staplers with Universal Handles			
Product	Code	Shaft Length	Pack QTY
Ultimate Endoscopic Linear Cutter Stapler with Universal Handle	PS9ECHU06	60mm	3 / Box
Ultimate Endoscopic Linear Cutter Stapler with Universal Handle	PS9ECHU16	160mm	3 / Box
Ultimate Endoscopic Linear Cutter Stapler with Universal Handle	PS9ECHU26	260mm	3 / Box

Ultimate Endoscopic Linear Cutter Stapler Articulating Reloading Units									
Product	Code	Colour	Tissue Type	Cut Length	Rows	Staples	Open Staple Height	Closed Staple Height	Pack QTY
Ultimate 30mm Endoscopic Articulating Reloading Unit	PS9ECAR30W	White	Thin	25mm	6	48	2.5mm	1.0mm	6 / Box
Ultimate 30mm Endoscopic Articulating Reloading Unit	PS9ECAR30B	Blue	Standard	25mm	6	48	3.5mm	1.5mm	6 / Box
Ultimate 45mm Endoscopic Articulating Reloading Unit	PS9ECAR45W	White	Thin	40mm	6	66	2.5mm	1.0mm	6 / Box
Ultimate 45mm Endoscopic Articulating Reloading Unit	PS9ECAR45B	Blue	Standard	40mm	6	66	3.5mm	1.5mm	6 / Box
Ultimate 45mm Endoscopic Articulating Reloading Unit	PS9ECAR45D	Gold	Regular - Thick	40mm	6	66	4.2mm	1.8mm	6 / Box
Ultimate 45mm Endoscopic Articulating Reloading Unit	PS9ECAR45G	Green	Thick	40mm	6	66	4.8mm	2.0mm	6 / Box
Ultimate 45mm Endoscopic Articulating Reloading Unit	PS9ECAR45N	Black	Very Thick	40mm	6	66	5.0mm	2.35mm	6 / Box
Ultimate 60mm Endoscopic Articulating Reloading Unit	PS9ECAR60W	White	Thin	55mm	6	90	2.5mm	1.0mm	6 / Box
Ultimate 60mm Endoscopic Articulating Reloading Unit	PS9ECAR60B	Blue	Standard	55mm	6	90	3.5mm	1.5mm	6 / Box
Ultimate 60mm Endoscopic Articulating Reloading Unit	PS9ECAR60D	Gold	Regular - Thick	55mm	6	90	4.2mm	1.8mm	6 / Box
Ultimate 60mm Endoscopic Articulating Reloading Unit	PS9ECAR60G	Green	Thick	55mm	6	90	4.8mm	2.0mm	6 / Box
Ultimate 60mm Endoscopic Articulating Reloading Unit	PS9ECAR60N	Black	Very Thick	55mm	6	90	5.0mm	2.35mm	6 / Box

Purple Surgical is a leading independent British manufacturer of high quality, cost-effective single-use surgical instrumentation and devices, partnering with global healthcare providers for positive patient outcomes.

Ultimate Endoscopic Linear Cutter Staplers with Universal Handles

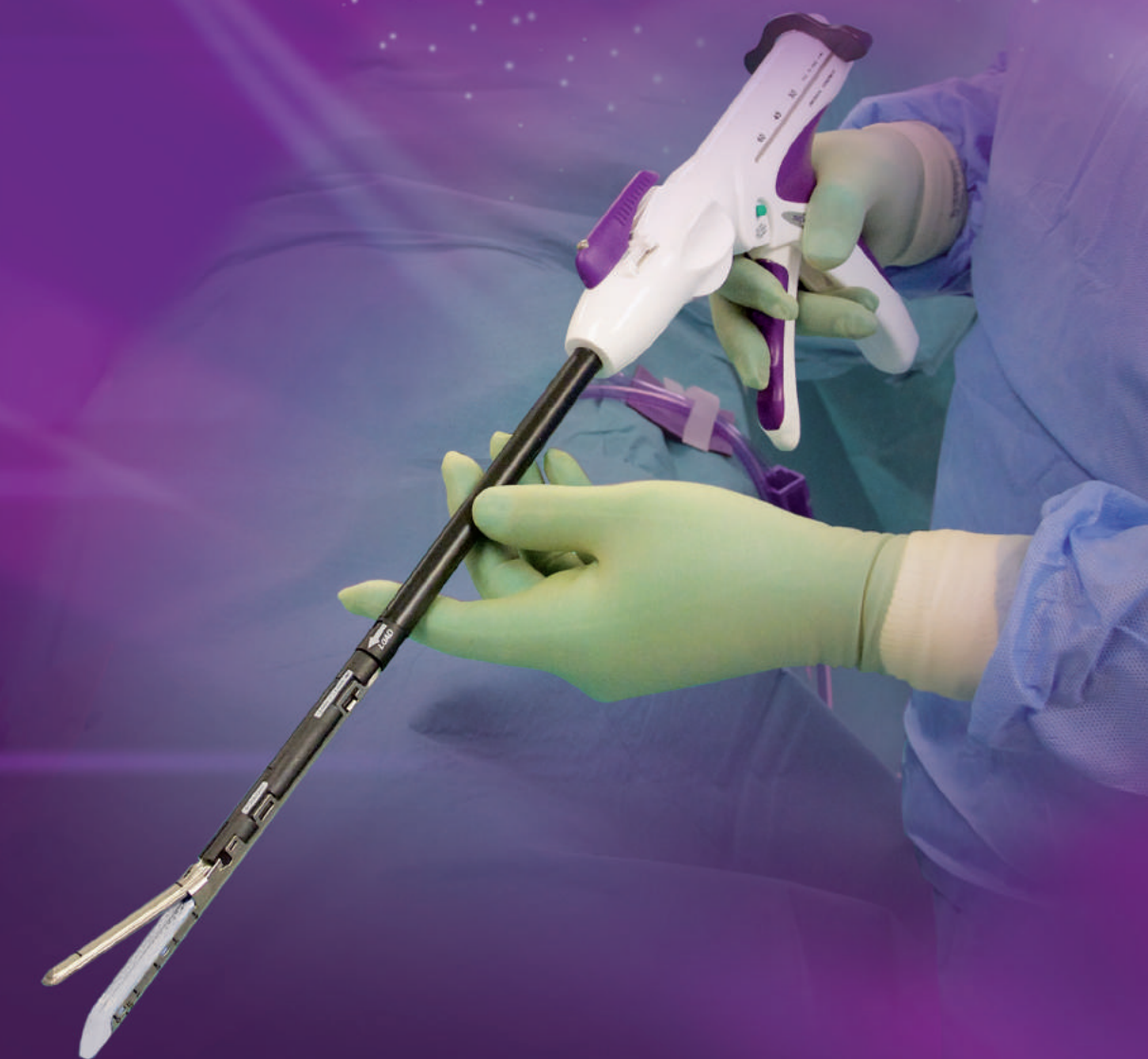

Contact Us
 Tel: +44 (0)1923 839333
 Email: mail@pursesurgical.com
www.pursesurgical.com



Great British Manufacturing
 Purple Surgical UK Ltd
 Purple Surgical International Ltd
 2 Chestnut House, Farm Close, Shenley, Herts. WD7 9AD, England

PS-ULT-ESTAP-22-09-17




FROM THE UK'S LEADING MANUFACTURER OF LAPAROSCOPIC INSTRUMENTS AND DEVICES

The Purple Surgical Ultimate Endoscopic Linear Cutter Staplers

The Purple Surgical *Ultimate* Endoscopic Linear Cutter Stapler places two, triple-staggered rows of titanium staples, and simultaneously cuts and divides the tissue between the two rows.

- The Universal Handle is available in 3 shaft lengths of 6cm, 16cm and 26cm, accommodating 30mm, 45mm and 60mm articulating reloading units.
- The reloading units are available in the following sizes, for varying tissue thicknesses:
 - White (2.5mm)
 - Blue (3.5mm)
 - Gold (4.2mm)
 - Green (4.8mm)
 - Black (5.0mm)
- White, Blue and Gold reloading units are compatible with Purple Surgical's *Ultimate* 12mm Trocars.
- Green and Black reloading units fit through 15mm Trocars.
- Each reloading unit houses a sharp knife blade to ensure a clean cut on every firing.
- Staple gap control allows for consistent staple formation as the cutting knife proceeds through the transecting tissue.

Articulating Reloading Units

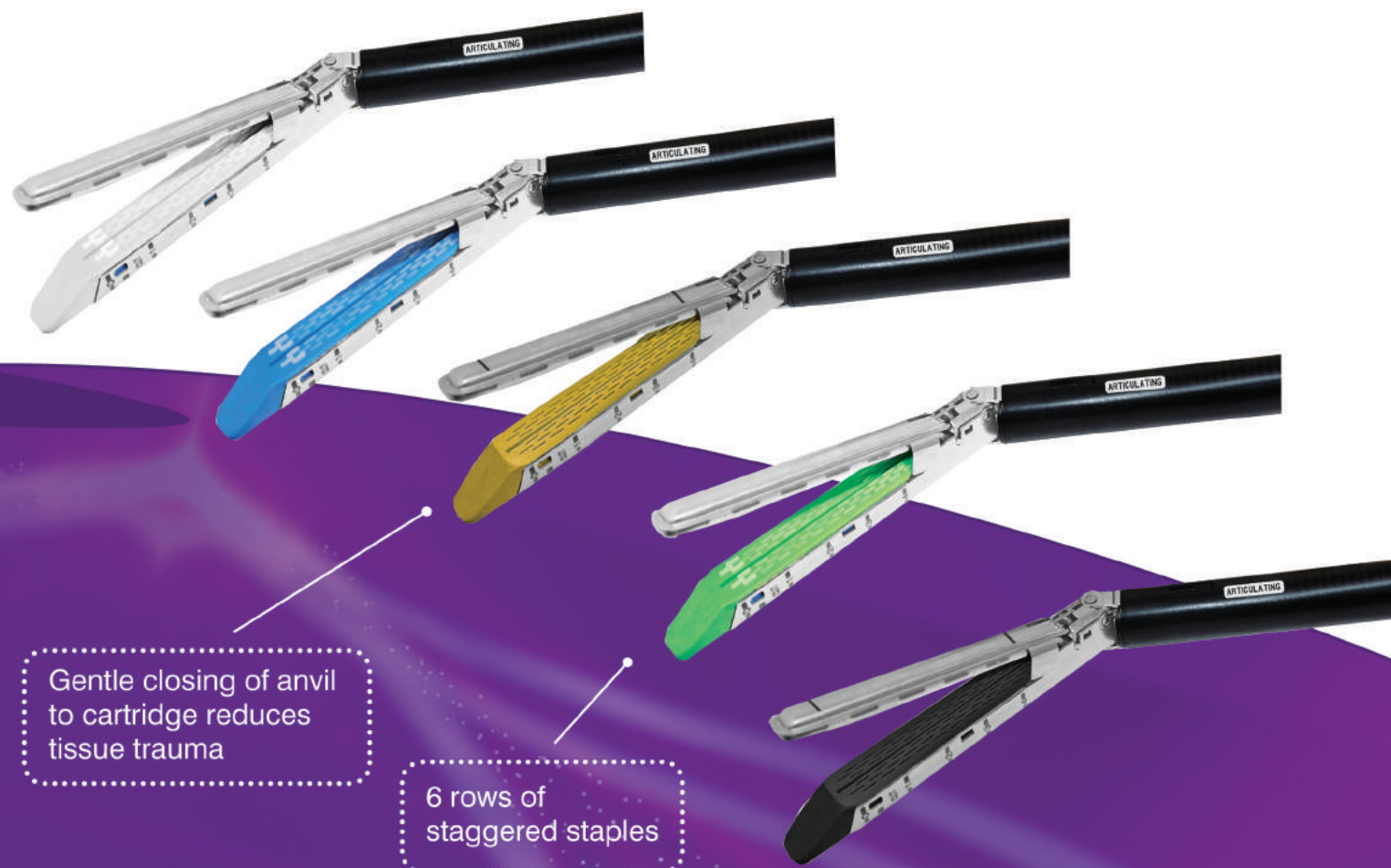


Gentle closing of anvil to cartridge reduces tissue trauma

6 rows of staggered staples

Features & Benefits:

- Gentle tissue compression before firing supports repositioning and minimises tissue trauma.
- The cartridge lock-out safety mechanism prevents the reuse of fired reloading units.
- The Universal Handle is available in three size options and is compatible with all Purple Surgical reloading units.
- The Universal Handle can be fired up to 25 times during a single procedure.
- Suitable for:
 - Bariatric Surgery
 - Colorectal Surgery
 - Gynaecological Surgery
 - Paediatric Surgery
 - Thoracic Surgery

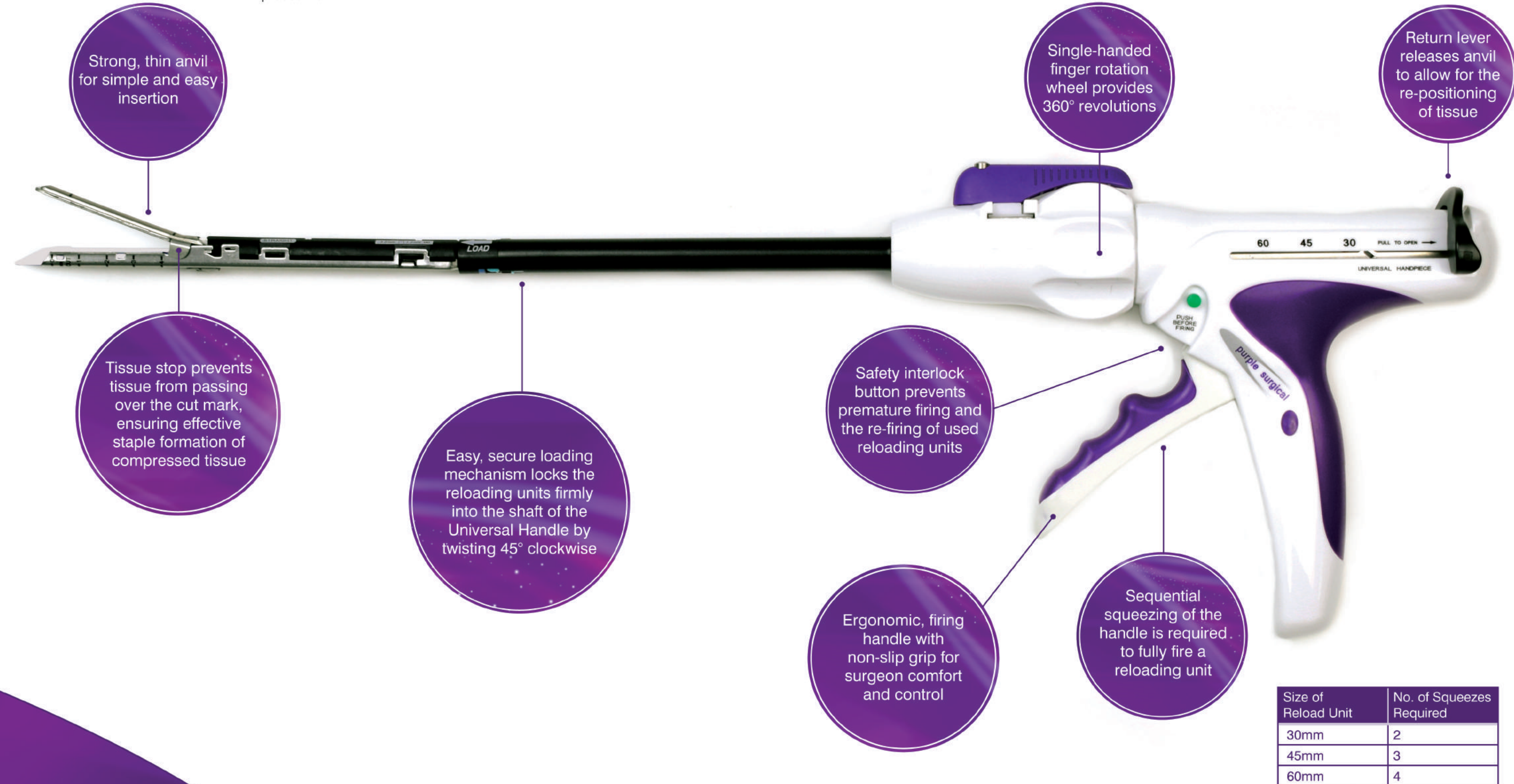

Strong, thin anvil for simple and easy insertion

Tissue stop prevents tissue from passing over the cut mark, ensuring effective staple formation of compressed tissue

Easy, secure loading mechanism locks the reloading units firmly into the shaft of the Universal Handle by twisting 45° clockwise

Single-handed finger rotation wheel provides 360° revolutions

Return lever releases anvil to allow for the re-positioning of tissue

Safety interlock button prevents premature firing and the re-firing of used reloading units

Ergonomic, firing handle with non-slip grip for surgeon comfort and control

Sequential squeezing of the handle is required to fully fire a reloading unit

Size of Reload Unit	No. of Squeezes Required
30mm	2
45mm	3
60mm	4

"Purple Surgical are hugely creative in finding the right solution for the challenges that hospitals are up against"

Liz Roberts, Theatre Stores Manager